

DTC3B N95

Product Technical Bulletin



Product Information

Model	Standard	Description	Color	Packing
DTC3B	NIOSH N95 FDA	Particulate Mask For Medical Use	White Fabric White headstrap	20 pcs/box 20 boxes/ctn

Product Features

- ✧ N95, Filters at least 95% of airborne particles. Not resistant to oil.
- ✧ Headstrap design, fit closer to the face.
- ✧ Adjustable nose clip for better fit.
- ✧ Ultrasonic welding technology, glue-free and odorless.

Product Materials

- ✧ Latex free
- ✧ Static-electro treated meltblown filtration material
- ✧ Soft woven elastic straps
- ✧ Adjustable Plastic Nosepiece and micro-aperture foam
- ✧ Non-toxic and odorless

Standard Information

Property	Typical Values
NIOSH	
N95 certification	TC-84A-4336
Filtration at 85 l/min	>95% of 0.3 micron particles
FDA	
510(K) NUMBER	K090131
Regulatory Class	II
Product Code	MSH, FXX
Listing Number	D198031

Container Loading QTY

✧ MEAS: 52 x 42 x 32 cm (0.069888cbm)

20' GP /28 cbm	40' GP/ 58 cbm	40' HQ / 68 cbm	45'HQ / 72cbm
400ctns (160,000 pcs)	825ctns (330,000 pcs)	970ctns (388,000 pcs)	1030ctns (412,000 pcs)

Package Pictures

Box inside




Insert certificate of conformity

合格证

产品名称: 口罩 (医用外科 N95)
 产品型号: DTC3B N95
 产品规格: 175×95mm
 主要成分: 无纺布 70%、熔喷布 30%
 执行标准: NIOSH N95
 过滤效率: ≥95%
 生产日期: 见盒底
 生产批号: 见盒底
 有效期: 3 年
 检验员: 01
 生产企业: 上海大胜卫生用品制造有限公司
 生产地址: 上海市松江区施惠路 228 号
 电话: 021-57783126 传真: 021-57784148

Box side

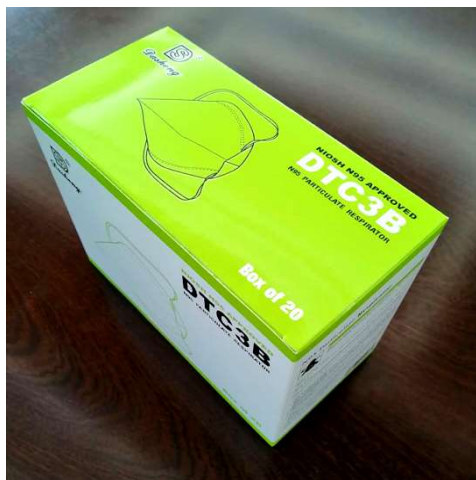


Box side



Package Artworks


✧ Box artwork

MEAS: 20 x 10x 15 cm

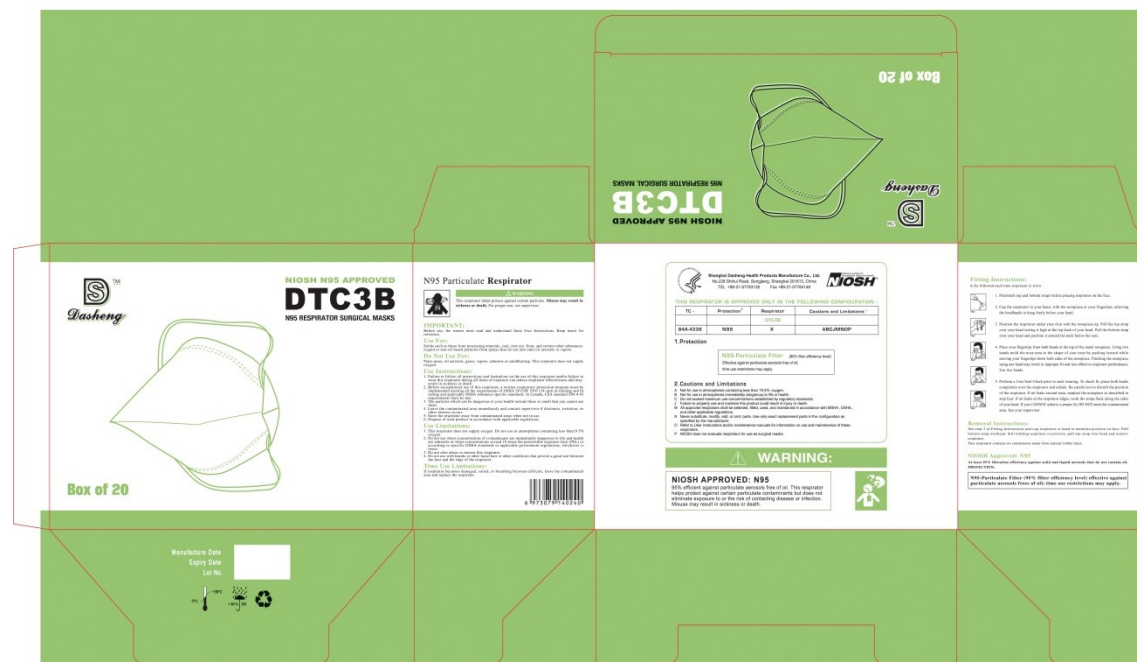

✧ Carton artwork

MEAS: 52 x 42x 32 cm

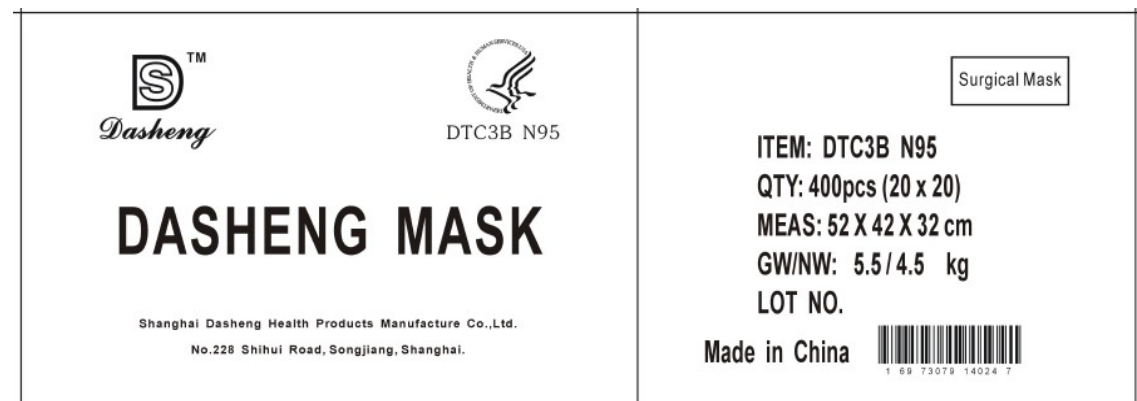



SHANGHAI DASHENGHEALTH
PRODUCTS MANUFACTURE CO LTD

Certificate of Conformity

✧ Insert each box

✧ Shall update time to time without notice

CERTIFICATE OF CONFORMITY

合格证

产品名称: 口罩 (医用外科 N95)

产品型号: DTC3B N95

产品规格: 175×95mm

主要成分: 无纺布 70%、熔喷布 30%

执行标准: NIOSH N95

过滤效率: ≥95%

生产日期: 见盒底

生产批号: 见盒底

有效期: 3 年

检验员: 01

生产企业: 上海大胜卫生用品制造有限公司

生产地址: 上海市松江区施惠路 228 号

电话: 021-57783126 传真: 021-57784148

Product Name: Surgical N95 Facemask

Model No.: DTC3B N95

Spec.: 175×95 mm

Main Material: Nonwoven, Meltblown

Executive Standard: NIOSH N95

Filtering Efficiency: ≥95%

Manufacture Date: see the bottom of the box

Lot No.: see the bottom of the box

Expiration Date: 3 years

Inspector: Stamp of QC PASS / No.

Manufacturer: Shanghai Dasheng Health Products
Manufacture Co., Ltd.

Address: No. 228 Shihui Road, Songjiang, Shanghai

Tel: 021-57783126 **Fax:** 021-57784148



SHANGHAI DASHENGHEALTH PRODUCTS MANUFACTURE CO LTD

NIOSH APPROVAL / DTC3B N95 TC-84A-4336

✧ N95 Approval Letter



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

NIOSH Reference: TN-14519
Mfr. Reference: SDHDT3B-8

Centers for Disease Control
and Prevention (CDC)
National Institute for Occupational
Safety and Health (NIOSH)
National Personal Protective
Technology Laboratory (NPPTL)
P.O. Box 18070
Pittsburgh, PA 15236-0070
Phone: 412-386-4000
Fax: 412-386-4051
August 2, 2006

Ms. Maggie Zhong
Shanghai Dasheng Health Products Manufacture Co., Ltd.
Room 604, No. 7 Building
No.20 Handan Road
Shanghai, 200437
CHINA

Dear Ms. Zhong:

The National Institute for Occupational Safety and Health (NIOSH) has reviewed your request dated May 5, 2006. This request was for an approval of the model DTC3B N95 filtering facepiece respirator.

This request is granted. Approvals are granted only for documentation written in the English language. It is the manufacturer's responsibility to correctly translate materials desired in languages other than English. Approval number TC-84A-4336 has been assigned. The respirator is approved for protection at an N95 particulate efficiency level.

The diskette enclosed with this letter contains the final respirator approval label. The abbreviated label has been accepted as submitted. The cautions and limitations, which apply to this approval, are on the approval label. Only those assemblies affected by this request, or where new approval numbers are assigned, apply to this approval action. Production approval labels cannot include information on unapproved configurations.

The approved assembly consists of the parts as listed on the approval labels plus the assembly matrices. Parts are to be marked with the numbers indicated on the approval label in a legible and permanent manner (marking cannot be removed without evidence of its previous presence). This certificate of approval is not an endorsement of the respirator by NIOSH, and such endorsement shall not be stated or implied in advertisements or other publicity. However, you may publicize the fact that this respirator has met the requirements of Title 42, *Code of Federal Regulations*, Part 84.



SHANGHAI DASHENGHEALTH
PRODUCTS MANUFACTURE CO LTD

Page 2 – Ms. Maggie Zhong – TN-14519

No additional changes may be made to any respirators and accompanying documentation without prior written approval of NIOSH. Requests for changes must be submitted to NIOSH and a modification of this approval must be granted before changes are made.

Sincerely yours,

Heinz W. Ahlers
Chief, Technology Evaluation Branch
National Personal Protective Technology Laboratory

Enclosures

✧ N95 Approval Label

On each box



Shanghai Dasheng Health Products Manufacture Co., Ltd.
No.228 Shihui Road, Songjiang, Shanghai 201613, China
TEL +86-21-57783126 Fax +86-21-57784148



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION :

TC -	Protection ¹	Respirator	Cautions and Limitations ²
		DTC3B	
84A-4336	N95	X	ABCJMNOP

1. Protection

N95-Particulate Filter (95% filter efficiency level)

Effective against particulate aerosols free of oil;
time use restrictions may apply.

2. Cautions and Limitations

- A Not for use in atmospheres containing less than 19.5% oxygen.
- B Not for use in atmospheres immediately dangerous to life or health.
- C Do not exceed maximum use concentrations established by regulatory standards.
- J Failure to properly use and maintain this product could result in injury or death.
- M All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O Refer to User Instructions and/or maintenance manuals for information on use and maintenance of these respirators.
- P NIOSH does not evaluate respirators for use as surgical masks.